

suPARnostic® TurbiLatex Reagents

Instructions for Use

REF T005

suPARnostic® and the ViroGates logo are registered trademarks of ViroGates A/S Denmark.
©2008 ViroGates.
All rights reserved.



CE IVD

ViroGates A/S
Banevænget 13
Birkerød 3460,
Denmark
Tel: +45 2113 1336
www.virogates.com

This product is protected by one or more US, European, and/or foreign patents.

The T005 product is validated on the automated Roche cobas c 701/2 instrument (cobas® is a trademark of Roche), and this instruction for use is dedicated to this biochemistry analyser.

Refer to the webpage <http://www.virogates.com> for instructions for use for other biochemistry analysers and other languages. Alternatively, contact your local distributor for instructions in your language.

INTENDED PURPOSE

For in vitro diagnostic use.

The suPARnostic® TurbiLatex Reagents are an in vitro diagnostic assay used to determine the soluble urokinase Plasminogen Activator Receptor (suPAR) level in human K2-EDTA and Lithium Heparin plasma on automated biochemistry analysers. The suPARnostic® TurbiLatex is a quantitative test measuring the suPAR level in ng/mL. It is intended to aid in the detection and evaluation of inflammatory disorders and immune activation.

INTENDED USER AND PATIENT

For professional use.

Typical users are laboratory technicians in central laboratories.

The typical patients are present at the Emergency Departments (ED) or the Intensive Care Units (ICU).

Acute Medicine

For unselected acute care patients, suPARnostic® TurbiLatex is used to identify the level of inflammation and immune activation to support triage decisions in conjunction with clinical findings and the results of other laboratory tests.

COVID-19

In patients with confirmed COVID-19 virus, suPARnostic® TurbiLatex is used to identify the level of inflammation and immune activation to aid in determining the risk of respiratory failure with a need for mechanical ventilation in conjunction with clinical findings and the results of other laboratory tests.

suPAR IS A MARKER OF DISEASE PROGRESSION

suPAR is the soluble form of urokinase Plasminogen Activator Receptor (uPAR). The amount of suPAR is a measure of immune activation and inflammation.¹ suPAR is a biomarker increased by disease presence and severity.

In unselected ED patients, suPAR has a high negative predictive value for ruling out disease progression². This means that patients with a low (<4 ng/mL) suPAR level have a good prognosis and a low risk of readmission and mortality³, supporting the decision of discharge of the patient. Conversely, a high suPAR level (>6 ng/ml) is a strong measure of chronic inflammation and the underlying risk of adverse outcomes, including short term mortality (in hospital, 30 days, or 90 days)² supporting the decision of further examination of the patient.

The use of suPAR in clinical routine adds significant additional knowledge to the standard assessment based on early warning score systems and standard parameters in pre-admission of acute medical patients. Hence, suPAR is a broadly applicable biomarker, e.g., in the ED, especially concerning discharge decisions of patients and identifying non-diagnosed inflammatory disease. A cluster-randomised interventional study showed that up- or down-triaging patients based on suPAR levels increased the number of patients discharged (low risk) by 34%⁴ and reduced hospital bed days⁵.

In patients with confirmed COVID-19, suPAR levels below 4 ng/mL suggest a low risk of developing respiratory failure and can be discharged for home quarantine.⁶

TEST PRINCIPLES

The suPARnostic® TurbiLatex test is a Particle Enhanced Turbidimetric Immunoassay (PETIA) that quantitatively determines suPAR in human plasma samples. The latex enhanced particles coated with anti-suPAR antibodies (Mouse/Rat) in the reagent agglutinates with the suPAR present in the sample. During the incubation time, an antigen-antibody complex is formed. The size of the complex is estimated using spectrophotometric technology at a wavelength of 570-590 nm range. The degree of turbidity caused by the agglutination is a measure of suPAR in the sample. The higher suPAR, the higher the turbidity.

REAGENTS AND MATERIALS

Reagents provided:

- Reagent 1: Dilution Buffer (Glycine-Buffer solution (pH 8.2) and preservatives)
- Reagent 2: Latex Particle Reagent (Phosphate buffer solution (pH 6.1), latex particles coated with anti-suPAR antibodies and preservatives)

This kit consists of a ready-to-use Reagent 1 dilution buffer and a ready-to-use Reagent 2 solution of latex particles coated with anti-suPAR antibodies.

The volume below is enough for a full cassette for the Roche cobas c 701 and c 702 chemistry analyser.

Product/Buffer	Reagent 1	Reagent 2	No. tests
T005 (c 701 or c 702)	75 mL	28 mL	420*

*No. of tests refers to the measurements available for sample testing when all reagent is transferred to the empty cassette. Reagent volumes additionally include the cassettes dead volumes and 12 tests used for calibration. The cassettes are for one-time use.

Material required but not provided:

- suPARnostic® TurbiLatex Calibrators
- suPARnostic® TurbiLatex Controls
- Roche Diagnostic cobas c 701 or c 702 analyser
- Empty cassettes for Roche cobas c 701 or c 702
- General laboratory equipment

PRECAUTIONS

For professional laboratory use.

For in vitro diagnostic use. Exercise the standard precautions required for handling all laboratory reagents. Disposal of all waste material should follow local guidelines. The safety data sheet is available for professional users upon request.

- Do not use kit components beyond the indicated kit expiration date.
- Do not freeze any of the kit components.
- Do not mix reagents from different kit lots.
- Do not switch caps on reagent containers as it may cause contamination or mix-up.
- Do not mouth pipette or ingest any of the reagents.
- Do not smoke, eat, or drink when performing the assay or in areas where samples or reagents are handled.
- Do not mix plasma samples from different patients or different blood samplings of the same patient.
- Human samples may be contaminated with infectious agents. Therefore, do not ingest, expose to open wounds, or breathe aerosols.
- Wear protective gloves and dispose of biological samples according to the regulations.
- Be aware of possible dilution of suPAR in the case of a transfusion, infusion or similar.

STORAGE AND HANDLING

The suPARnostic® TurbiLatex Reagents kit should be stored at 2–8°C – do not freeze.

Before use, check the expiry date on the label.

The suPARnostic® Reagents are produced with a 2-year shelf life from the production date.

The reagents have 8 weeks of onboard stability when kept at 2–8°C and with a minimum of monthly calibration.

The reagent's stability may be affected if not stored correctly, and inefficient and misleading results may be obtained. If any colouring or precipitation appear, discard the Reagents.

The suPARnostic® TurbiLatex Reagent must be transferred to Roche's appropriate cassettes before loading on the equipment.

SAMPLE COLLECTION AND PREPARATION

Blood samples collected in K2-EDTA and Lithium-Heparin have been validated.

Collecting blood samples should be performed by authorised personnel using approved venipuncture techniques.

To prepare plasma samples, draw whole blood into a blood collection tube containing either K2-EDTA or Lithium Heparin anticoagulant. Then, centrifuge the blood at 3,000 x g between 1–10 minutes or until blood cells and plasma have separated.

Ensure the samples, calibrators, and controls are at room temperature before measurement.

Due to possible evaporation effects when loaded on the instrument, samples, calibrators, and controls should be analysed within 2 hours.

NOTE: Do not use hemolysed, contaminated, or hyper-lipemic sample specimens.

ASSAY PROCEDURE

- 1) Install the suPARnostic® method using the application parameters provided at the end of this note onto the Roche cobas c 701/2 analyser.
- 2) Transfer the suPARnostic® TurbiLatex Reagents to Roche's appropriate cassettes before loading on the equipment. Reagent 1 is transferred to position A and Reagent 2 to position C. For more information, refer to the Roche cobas c MULTI pack IFU.

NOTE: The suPARnostic® TurbiLatex Reagent 2 is on the Roche cobas platform identified as R3.

- 3) Load reagents on the analyser.
- 4) Load blood sample directly on the analyser if fully automated or isolate plasma before loading.
- 5) The assay run time is 10 minutes with:
 - 1st incubation: 150 µL of Reagent 1 with 10 µL of the sample.
 - 2nd incubation: 50 µL of Reagent 2 is added to the mixture, and antigen-antibody complexes are formed.
- 6) The turbidity of the sample is measured in set time intervals at a wavelength of 570 nm.
- 7) Results are determined via a calibration curve generated by measuring a set of calibrators (#T002) with a known suPAR concentration.
- 8) The measurement result is calculated by determining the difference of absorbance values at 2 read points. The analyser automatically calculates the analyte concentration of each sample in ng/mL.

CALIBRATION

Together with the suPARnostic® TurbiLatex Reagents kit, the suPARnostic® TurbiLatex Calibrators (#T002) must be used for calibration. It is recommended to repeat the calibration at least once a month. In addition, recalibration is required when a new batch of the suPARnostic® TurbiLatex Reagents are used. Calibration is performed according to the instruction provided with the suPARnostic® TurbiLatex Calibrators.

QUALITY CONTROL

Quality control of the suPARnostic® TurbiLatex Reagents must be performed with the suPARnostic® TurbiLatex Controls (#T003) as a minimum after each calibration and according to the laboratory guidelines. If the QC values exceed the established upper/lower range, the laboratory should undertake corrective actions.

INTERPRETATION OF RESULTS

suPAR levels and cut-offs

Acute care medical patients and risk of 90-day mortality

The cut-offs for interpreting results from acute care patients were established based on suPAR baseline measurements of 990 patients attending the ED in a Spanish multicenter trial.¹⁴ The median age was 68 years (53-81), 50.8% were men, median suPAR was 3.8 ng/mL (interquartile range 2.8-6.0). In total, 47 died during the 90-day follow-up. Out of the 990 patients, 520 (52.5%)

had suPAR below 4.0 ng/mL. Patients with suPAR <4.0 ng/mL had a low risk of 90-day mortality (N=5, 0.96%), resulting in a negative predictive value (NPV) of 99.0%, a sensitivity of 89.4% and a specificity of 54.6%. In patients with suPAR >6.0 ng/mL (N=245 (24.8%)), 33 patients died during the 90-day follow-up (13.5%), resulting in a positive predictive value (PPV) of 13.5%, a sensitivity of 70.2%, and a specificity of 77.5%.

	90 Days follow-up		Total	PPV	NPV
	Died	Survived			
High risk (suPAR >6.0 ng/mL)	33	212	245	13.5%	
Medium Risk (suPAR 4.0– 6.0 ng/mL)	9	216	225		
Low risk (suPAR <4.0 ng/mL)	5	515	520		99.0%
Total	47	943	990		
Sensitivity/specificity (<4.0 ng/mL)	89.4%	54.6%			
Sensitivity/specificity (>6.0 ng/mL)	70.2%	77.5%			

Table 1: 90-day Mortality according to suPAR cut-offs in Spanish multicenter study.

COVID-19 and risk of respiratory failure

For patients who tested positive for the COVID-19 virus, suPAR baseline measurements were taken within 48 hours after patients presented at the hospital⁶. Respiratory failure was defined as the need for mechanical ventilation within 2 weeks. The study included 57 patients, of whom 21 developed respiratory failure. None of the patients with suPAR below 4.0 ng/mL developed respiratory failure resulting in an NPV of 100%, a sensitivity of 100% and a specificity of 36.1%. Of the 21 patients who developed respiratory failure, 18 had baseline suPAR levels above 6.0 ng/mL resulting in a PPV of 85.7%, a sensitivity of 85.7%, and a specificity of 81.3%.

suPAR level	Interpretation, ED and COVID-19
<4.0 ng/mL	<p>Low Risk</p> <ul style="list-style-type: none"> - Supports the decision of discharge. - The underlying health condition is good, and the prognosis for survival is high. - Low risk of respiratory failure and need of mechanical ventilation in patients with COVID-19.
4.0-6.0 ng/mL	<p>Medium Risk</p> <ul style="list-style-type: none"> - Some disease activity or co-morbidity is present. - Some readmissions and mortality are expected after six months of follow-up. - Medium risk of respiratory failure and need of mechanical ventilation in patients with COVID-19.
>6.0 ng/mL	<p>High Risk</p> <ul style="list-style-type: none"> - Clinical attention is needed - high risk of mortality. - Supports the decision of admission and treatment - High risk of respiratory failure and need of mechanical ventilation in patients with COVID-19.

Table 2: Simplified suPAR clinical decision scheme^{6,14}.

EXPECTED VALUES IN HEALTHY INDIVIDUALS

All individuals have a measurable suPAR level. In healthy blood donors (N=9305), the median suPAR level for men aged 18-65 years is 2.2 ng/mL (25-75% interval from 1.8-2.9 ng/mL)⁷, for women aged

18–65 years is 2.6 ng/mL (25–75% interval from 2.1–3.2 ng/mL)⁷, patients attending ED the suPAR level is around 3.0–6.0 ng/mL^{2,3,8}. In patients with severe disease and organ failure, suPAR is often in the double-digits^{9,10}. The higher the level, the higher the risk of disease progression and the worse the prognosis.

CLINICAL PERFORMANCE

Validation of cut-offs

Acute care medical patients

The clinical validation data are from a prospective observational study of unselected acute medical patients attending the ED at Mikkeli Hospital in Finland.¹¹ A total of 1747 acute medical patients were included and had suPAR measured using suPARnostic® TurbiLatex. The median age was 70 (IQR: 57–79), and 51.4% were men. Among patients with suPAR below 4.0 ng/mL (N=804, 46.0%), 8 (1.0%) died within 90-day follow-up, resulting in a negative predictive value of 99.0% and a sensitivity of 94.2% and a specificity of 47.9%. Among patients with suPAR above 6.0 ng/mL (N=429, 24.6%), 87 patients (20.3%) died within 90-day follow-up, resulting in a positive predictive value of 20.1%, a sensitivity of 63.0% and a specificity of 78.7%. Data for 90-day follow-up are shown in Table 3.

	90-day follow-up		Total	PPV	NPV
	Died	Survived			
High risk (suPAR >6.0 ng/mL)	87	342	429	20.3%	
Medium Risk (suPAR 4.0– 6.0 ng/mL)	43	471	514		
Low risk (suPAR <4.0 ng/mL)	8	796	804		99.0%
Total	138	1609	1.747		
Sensitivity/specificity (< 4.0 ng/mL)	94.2%	49.5%			
Sensitivity/specificity (> 6.0 ng/mL)	63.0%	78.7%			

Table 3: 90-day mortality in acute medical patients in Finish validation study.

COVID-19

The clinical validation data are from a prospective observational study at Mikkeli Central Hospital in Finland using the suPARnostic® TurbiLatex on a cobas c 701. The study included 100 acute medical patients who tested positive for SARS-CoV-2 at the Emergency Department at Mikkeli Central Hospital in Finland.¹⁵

Results of suPAR validation for stratification of COVID-19 patients regarding the risk of developing severe respiratory failure and requiring mechanical ventilation is shown in Table 4.

	90-day follow-up		Total	PPV	NPV
	Died	Survived			
High risk (suPAR >6.0 ng/mL)	5	44	49	10.2%	
Medium Risk (suPAR 4.0– 6.0 ng/mL)	0	27	27		
Low risk (suPAR <4.0 ng/mL)	0	24	24		100%
Total	5	95	100		
Sensitivity/specificity (< 4.0 ng/mL)	100%	25.3%			
Sensitivity/specificity (> 6.0 ng/mL)	100%	53.7%			

Table 4: 90-day development of respiratory failure in COVID-19 patients according to suPAR cut-offs.

LIMITATIONS

Clinical prognostication must not be based on the results of the suPARnostic® TurbiLatex test alone. Instead, the results must be interpreted considering the patient's clinical history and other diagnostic test results.

ANALYTICAL PERFORMANCE**SAMPLE STABILITY**

Blood samples should be added to the instrument for automated sampling within 2 hours of sampling to avoid hemolysis.

Samples should preferably be analysed as soon as possible, but K2-EDTA and Lithium Heparin plasma samples are stable for:

- 24 hours at room temperature (20-25°C).
- 3 days at 2-8°C.
- -20°C for more extended storage.

Five freeze/thaw cycles over 5 days do not impact suPAR concentration in a sample.

REQUIRED TRAINING

To use the suPARnostic® TurbiLatex Reagents, the user must be fully trained to operate the chemistry analyser.

PERFORMANCE CHARACTERISTICS

Results presented below were obtained with the use of the suPARnostic® TurbiLatex Reagents on the Roche cobas® c 702 analyser. Therefore, the data shown is only valid for Roche cobas c 701 and c 702.

RESULTS

Results are calculated by linear regression.

If a calculation method other than linear regression is used, the method must be validated. Then, control the curve fitting using the suPARnostic® TurbiLatex Controls and undertake corrective actions if results exceed the upper and lower limit.

MEASURING RANGE

The measuring range of the suPARnostic® TurbiLatex test is 1.8 ng/mL to 16.0 ng/mL on the Roche cobas® c 701/c 702 analyser.

It is not recommended to dilute samples with the results above the measuring range.

ANALYTICAL LIMITS

The Limit of Blank (LoB) was determined as the 95th percentile of 60 blank measurements.

The Limit of Detection (LoD) was calculated as LoB plus 3 times the pooled standard deviation (SD).

The LoD corresponds to the lowest analyte concentration, which can be detected.

The Limit of Quantitation (LoQ) was calculated as LoB plus 10 times the SD of the mean low-level 'samples' estimation of 60 low-level samples.

LoB	LoD	LoQ
1.0 ng/mL	1.3 ng/mL	1.3 ng/mL

LoB and LoD were established following CLSI EPI7¹³. LoQ was established according to IUPAC, Codex Alimentarius Procedural Manual, 15th edition.

INTERFERENCE

Samples with abnormally elevated haemoglobin levels, lipids, or bilirubin may interfere with assay performance and sensitivity.

No interference was observed for the following concentrations.

Substance:	Concentration:
Bilirubin	350 µmol/L
Hemoglobin	1.4 g/L
Triglycerides	3.3 g/L
Rheumatoid Factor	>440 IU/mL
HAMA	Titer >640*

The interference studies were performed using a modified CLSI EP07-A2 protocol.¹³

Rheumatoid Factor and HAMA solutions were prepared by adding concentrated rheumatoid and HAMA solutions to human plasma pools.

In rare cases, gammopathy, especially type IgM (Waldenström's macroglobulinemia), may produce inaccurate results. Patients with a confirmed diagnosis of anti-TPO or other autoimmune-related diseases have been shown to interfere in a few cases.

Although precautions have been taken to minimise interference caused by heterophilic antibodies, erroneous results can be observed. Therefore, any suPAR value above 10 ng/mL should be carefully investigated, and unusually high results, e.g. above 20 ng/mL may be false-positive results caused by interference.

LINEARITY

The suPARnostic® TurbiLatex Reagents test is linear from 1.8 ng/mL to 26.5 ng/mL.

HOOK EFFECT

The suPARnostic® TurbiLatex Reagents test showed no prozone effect in concentrations up to 47.5 ng/mL.

PRECISION

Low, medium, and high samples were measured with two replicates in two separate runs per day for 20 days.

	Mean suPAR level (ng/mL)	Repeatability (CV)	Within-day precision (CV)	Between-day precision (CV)	Within-laboratory precision (CV)
Low	3.7	3.5%	1.0%	10.8%	11.4%
Medium	6.9	3.5%	2.4%	4.6%	5.8%
High	9.9	4.1%	1.6%	4.5%	6.1%

The intermediate precision study was performed according to the CLSI EP05-A2 protocol.¹³

ACCURACY (METHOD COMPARISON)

The suPARnostic® TurbiLatex is calibrated against an internal control verified with suPARnostic® ELISA. The maximum allowed variation between suPARnostic® ELISA and TurbiLatex is 15%, and between lots of TurbiLatex, it is 10%.

Bias and correlation calculations toward the suPARnostic® ELISA kit were conducted to estimate the suPARnostic® TurbiLatex kits ability to quantify suPAR in patient samples. 103 samples were measured with one lot of suPARnostic® TurbiLatex reagents, and the results were compared to the suPARnostic® ELISA results.

Results:

Sample Type	No. of pairs	Slope	y-Intercept	Pearson correl.	Range value
Plasma	103	1.0957	-0.1036	0.959	1.8-15.7 ng/mL

X = suPARnostic® ELISA

Y = suPARnostic® TurbiLatex

ANTICOAGULANT EFFECT (METHOD COMPARISON)

The suPARnostic® TurbiLatex is calibrated for plasma samples with K2-EDTA anticoagulant. Therefore, bias and correlation calculations toward K2-EDTA based plasma samples should be considered when a Lithium Heparin anticoagulant is used.

Therefore, 45 samples from the same individual were collected in K2-EDTA and Lithium Heparin samples drawn from a single subject were measured with one lot of suPARnostic® TurbiLatex reagents, and the results were compared.

Sample Type	No. of pairs	Slope	y-Intercept	Pearson correl.	Range value
Lithium Heparin Plasma	45	1.1219	-0.0571	0.982	2.2-12.1 ng/mL

X = suPARnostic® TurbiLatex on K2-EDTA plasma.



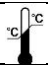
Y = suPARnostic® TurbiLatex on Lithium Heparin Plasma. Studies were performed using Siemens Atellica® bioanalyser.

When using plasma with Lithium-Heparin anticoagulant, the following correction factor must be added to offset the inherent matrix:

Correlated result = (Obtained result for Lithium Heparin plasma – 0.0571)/1.1219

WASTE HANDLING

Discard unused reagents and waste according to country, federal, state, and local regulations.

REF		
Catalogue no.	Contains sufficient for <n> tests	Use by
IVD		LOT
In vitro diagnostic medical device	Temperature Limits	LOT no. (Batch No.)

REFERENCES

- 1) Desmedt S et al. The intriguing role of soluble urokinase receptor in inflammatory diseases. *Crit Rev Clin Lab Sci.* 2017 Mar;54(2):117-133
- 2) Rasmussen LJH et al. Combining National Early Warning Score with Soluble Urokinase Plasminogen Activator Receptor (suPAR) Improves Risk Prediction in Acute Medical Patients: A Registry-Based Cohort Study. *Crit Care Med.* 2018 Dec;46:1961-8.
- 3) Rasmussen LJH, et al. Soluble urokinase plasminogen activator receptor (suPAR) in acute care: a strong marker of disease presence and severity, readmission and mortality. A retrospective cohort study. *Emerg Med J.* 2016 Nov;33:769-75.
- 4) Schultz et al. Availability of suPAR in emergency departments may improve risk stratification: a secondary analysis of the TRIAGE III trial *Scandinavian Journal of Trauma, Resuscitation and Emergency Medicine BMC* (2019) 27:43
- 5) Schultz, M. et al. Early Discharge from the Emergency Department Based on Soluble Urokinase Plasminogen Activator Receptor (suPAR) Levels: A TRIAGE III Substudy, Hindawi, *Disease Markers*, Volume 2019
- 6) Rovina et al. Soluble urokinase plasminogen activator receptor (suPAR) as an early predictor of severe respiratory failure in patients with COVID-19 pneumonia. *Crit Care*, 2020 Apr 30;24(1):187;
- 7) Haastrup E, et al.: Soluble urokinase plasminogen activator receptor as a marker for use of antidepressants. *PLoS One* 2014.
- 8) Raggam RB et al. Soluble urokinase plasminogen activator receptor predicts mortality in patients with systemic inflammatory response syndrome. *J Intern Med* 2014, 276(6):651-8
- 9) Koch A, et al. Clinical relevance and cellular source of elevated soluble urokinase plasminogen activator receptor (suPAR) in acute liver failure. *Liver Int* 2014;34:1330-1339.
- 10) Donadello K, et al. Soluble urokinase-type plasminogen activator receptor as a prognostic biomarker in critically ill patients. *J Crit Care.* 2014 Feb;29(1):144-9.
- 11) Seppälä, S. et al. suPAR Cut-offs for Stratification of Low, Medium, and High-risk Acute Medical Patients in the Emergency Department, preprint available at <https://www.researchsquare.com/article/rs-542503/v1>
- 12) Azam Tu, et al. International Study of Inflammation in COVID-19. Soluble Urokinase Receptor (SuPAR) in COVID-19-Related AKI. *J Am Soc Nephrol.* 2020 Nov;31(11):2725-2735.
- 13) Protocols for Determination of Limits of Detection and Limits of Quantitation; Approved Guideline. EPI7-A, Vol. 24 No. 34 Replaces EPI7-P, Vol. 24 No. 10. <https://www.clsi.org/>
- 14) Unpublished data from a multicenter study in Spain
- 15) Altintas I, et al.. suPAR Cut-Offs for Risk Stratification in Patients With Symptoms of COVID-19. *Biomark Insights.* 2021 Aug

APPLICATION PARAMETERS

For the Roche cobas c 701/2 instrument

Analyse		Calib.		Range			Other		
Assay/Time/Point	2Point End		10	22	37	0	0		
Wavelength (2./Pri.)	800		570						
Sample volume				Cassette configuration					
Norm.	10.0	0.0	0	code #####					
Dec.	10.0	0.0	0	Expiration Days 60					
inc.	10.0	0.0	0	Reagent volume					
Dilution				R1	150	0	Inactive		
<input type="checkbox"/> Water				R2	0	0	Inactive		
<input type="checkbox"/> Diluent				R3	50	0	Inactive		
Linearity Limit	0	%	0	%	0	0			
Prozone-Limit	-32000		32000	0	0	0	Inside	0	0
Abs. Limit	32000		Increase						
Cell Detergent	Detergent 1		Stirring Level			4			
Stirring Setting			M1	M2	M3				
UP	Stirring	LOW	Stirring	Stirring	Stirring				

Analyse	Calib.	Range	Other
Calibration type Linear		Auto-calibration	
Point	<input type="text" value="6"/>	<input type="checkbox"/> Timeout	
Span	<input type="text" value="3"/>	cassette	
Weight	<input type="text" value="0"/>	cassette	
Update Type	<input type="text" value="None"/> <input type="text" value="0"/> <input type="text" value="0"/>	<input type="text" value="Active"/>	
SD-Limit	<input type="text" value="999"/>	<input type="text" value="30"/> <input type="text" value="Day"/>	
Duplicate Limit	<input type="text" value="99"/> % <input type="text" value="32000"/> Abs.	Change over	
Sensitivity Limit	<input type="text" value="-99999"/> <input type="text" value="99999"/>	cassette <input type="text" value="Cancel"/>	
S1 Abs. Limit	<input type="text" value="-32000"/> <input type="text" value="32000"/>	<input type="checkbox"/> QC-Violation	
<input type="checkbox"/> Auto masking		Method <input type="text" value="Blank"/> Rule <input type="text" value="Is"/> Control1 <input type="text" value="None"/> Control2 <input type="text" value="None"/> Control3 <input type="text" value="None"/>	

Analyse	Calib.	Range	Other																																				
Applications code	###	<table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th colspan="4">Expected Values</th> </tr> <tr> <th colspan="4">Male</th> </tr> <tr> <td></td> <td></td> <td>-99999</td> <td>999999</td> </tr> <tr> <td>99</td> <td>Year</td> <td>-99999</td> <td>999999</td> </tr> <tr> <td>100</td> <td>Year</td> <td>-99999</td> <td>999999</td> </tr> <tr> <th colspan="4">Female</th> </tr> <tr> <td></td> <td></td> <td>-99999</td> <td>999999</td> </tr> <tr> <td>99</td> <td>Year</td> <td>-99999</td> <td>999999</td> </tr> <tr> <td>100</td> <td>Year</td> <td>-99999</td> <td>999999</td> </tr> </table>		Expected Values				Male						-99999	999999	99	Year	-99999	999999	100	Year	-99999	999999	Female						-99999	999999	99	Year	-99999	999999	100	Year	-99999	999999
Expected Values																																							
Male																																							
				-99999	999999																																		
99	Year			-99999	999999																																		
100	Year			-99999	999999																																		
Female																																							
				-99999	999999																																		
99	Year			-99999	999999																																		
100	Year			-99999	999999																																		
Unit	ng/mL																																						
Report Name	suPAR																																						
Data mode	Active																																						
<input type="checkbox"/> Automatic Rerun																																							
Technical Limit	1.8 16.0																																						
Repeat Limit	- 99999																																						
<input type="checkbox"/> Control interval Time	0																																						
<input type="checkbox"/> Automatic QC on-board- 1																																							
<input type="checkbox"/> Quality <table border="1" style="width:100%; border-collapse: collapse;"> <tr> <td>(1)</td> <td>0</td> <td>L</td> <td>0</td> </tr> <tr> <td>(2)</td> <td>0</td> <td>H</td> <td>0</td> </tr> <tr> <td>(3)</td> <td>0</td> <td>I</td> <td>0</td> </tr> <tr> <td>(4)</td> <td>0</td> <td></td> <td></td> </tr> <tr> <td>(5)</td> <td>0</td> <td></td> <td></td> </tr> <tr> <td>(6)</td> <td>0</td> <td></td> <td></td> </tr> </table>		(1)	0	L	0	(2)	0	H	0	(3)	0	I	0	(4)	0			(5)	0			(6)	0			<table border="1" style="width:100%; border-collapse: collapse;"> <tr> <th colspan="2">Default</th> </tr> <tr> <td colspan="2">Sex</td> </tr> <tr> <td><input type="checkbox"/> Male</td> <td><input type="checkbox"/> Female</td> </tr> <tr> <th colspan="2">Interval</th> </tr> <tr> <td><input type="checkbox"/> Range 1</td> <td><input type="checkbox"/> Range 2 <input type="checkbox"/> Range 3</td> </tr> </table>	Default		Sex		<input type="checkbox"/> Male	<input type="checkbox"/> Female	Interval		<input type="checkbox"/> Range 1	<input type="checkbox"/> Range 2 <input type="checkbox"/> Range 3			
(1)	0	L	0																																				
(2)	0	H	0																																				
(3)	0	I	0																																				
(4)	0																																						
(5)	0																																						
(6)	0																																						
Default																																							
Sex																																							
<input type="checkbox"/> Male	<input type="checkbox"/> Female																																						
Interval																																							
<input type="checkbox"/> Range 1	<input type="checkbox"/> Range 2 <input type="checkbox"/> Range 3																																						

Analyse	Calib.	Range	Other			
Standards	(1)	(2)	(3)	(4)	(5)	(6)
Calibrator code	*1	*1	*1	*1	*1	*1
Concentration	*2	*2	*2	*2	*2	*2
Rack-nr. pos.	*3	*3	*3	*3	*3	*3
Sample Volume	10.0	10.0	10.0	10.0	10.0	10.0
Diluted Sample volume	0.0	0.0	0.0	0.0	0.0	0.0
Diluent Volume	0	0	0	0	0	0

Cassette type for cobas® 6000 (c701)
Type A

Bottle			
a	N/A	N/A	N/A
b	R1	434	75
c	R3	434	28

Cassette type for cobas® 8000 (c702)
Type A

Bottle			
a	N/A	N/A	N/A
b	R1	434	75
c	R3	434	28

***1: User-defined value**

***2: Batch-specific concentrations - see CoA**

***3: User defined value**